

SynthMicro — Whole Blood Smear XAI Report


Analysis ID: 501ed74e239a4b8e83487169f5b0ada4
Generated: 2026-09-07 23:36:44

Pipeline: Cellpose segmentation → 384×384 cell extraction → ResNet50 classification → Grad-CAM → SHAP → explanatory report.

Measurement	Result
Cellpose objects	14
Nucleated-cell candidates	1
Myeloblast-like candidates	0
Myeloblast-like proportion	0.00%

Screening interpretation: the proportion above is a research-oriented model indicator based on candidate cells selected by the notebook's purple-score rule. It is not a cancer probability and does not establish AML or another diagnosis.


Whole Smear

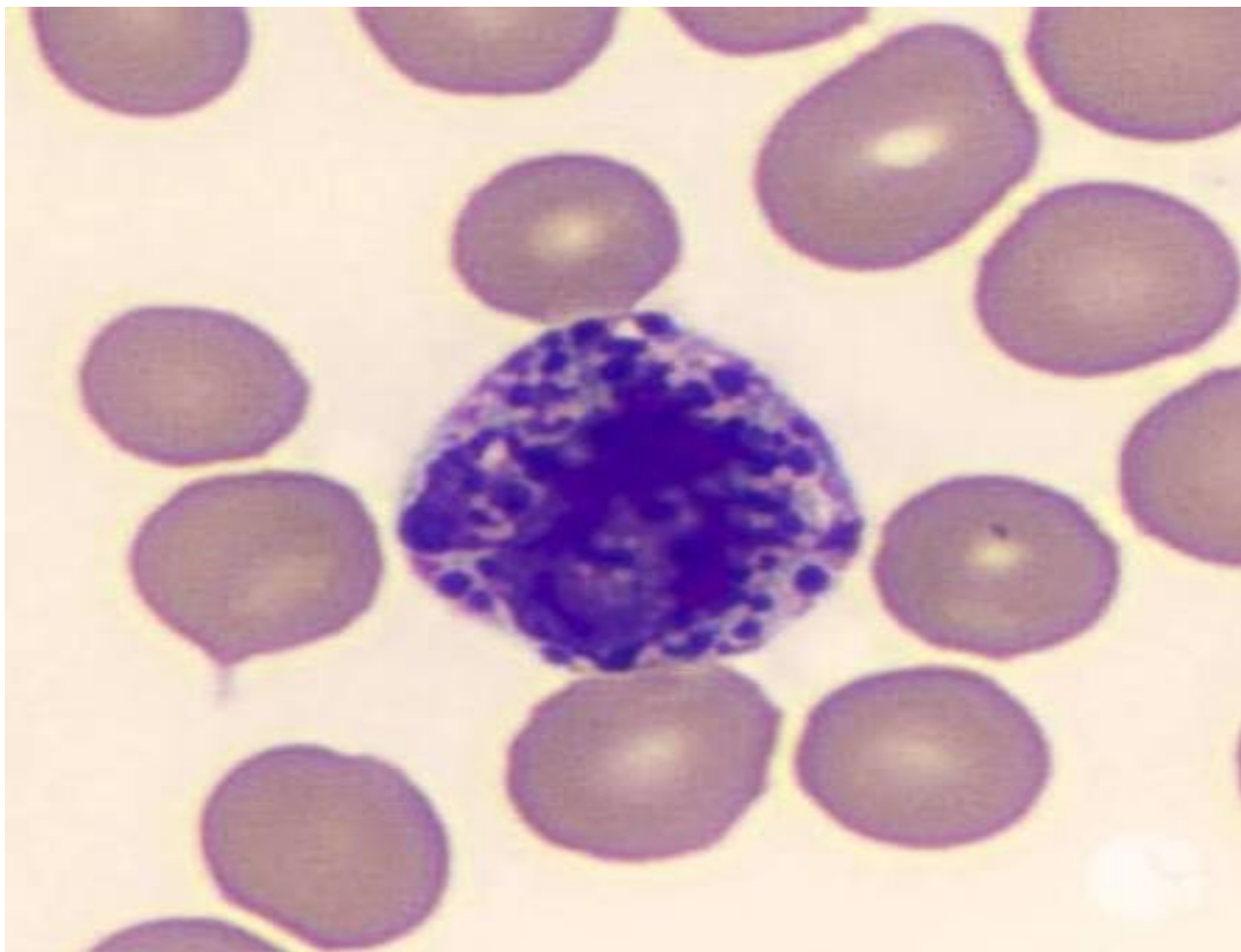

Original uploaded blood-smear image analyzed by Cellpose.

Segmentation and Candidate Overlay

Whole Smear — Candidate Cell Classification

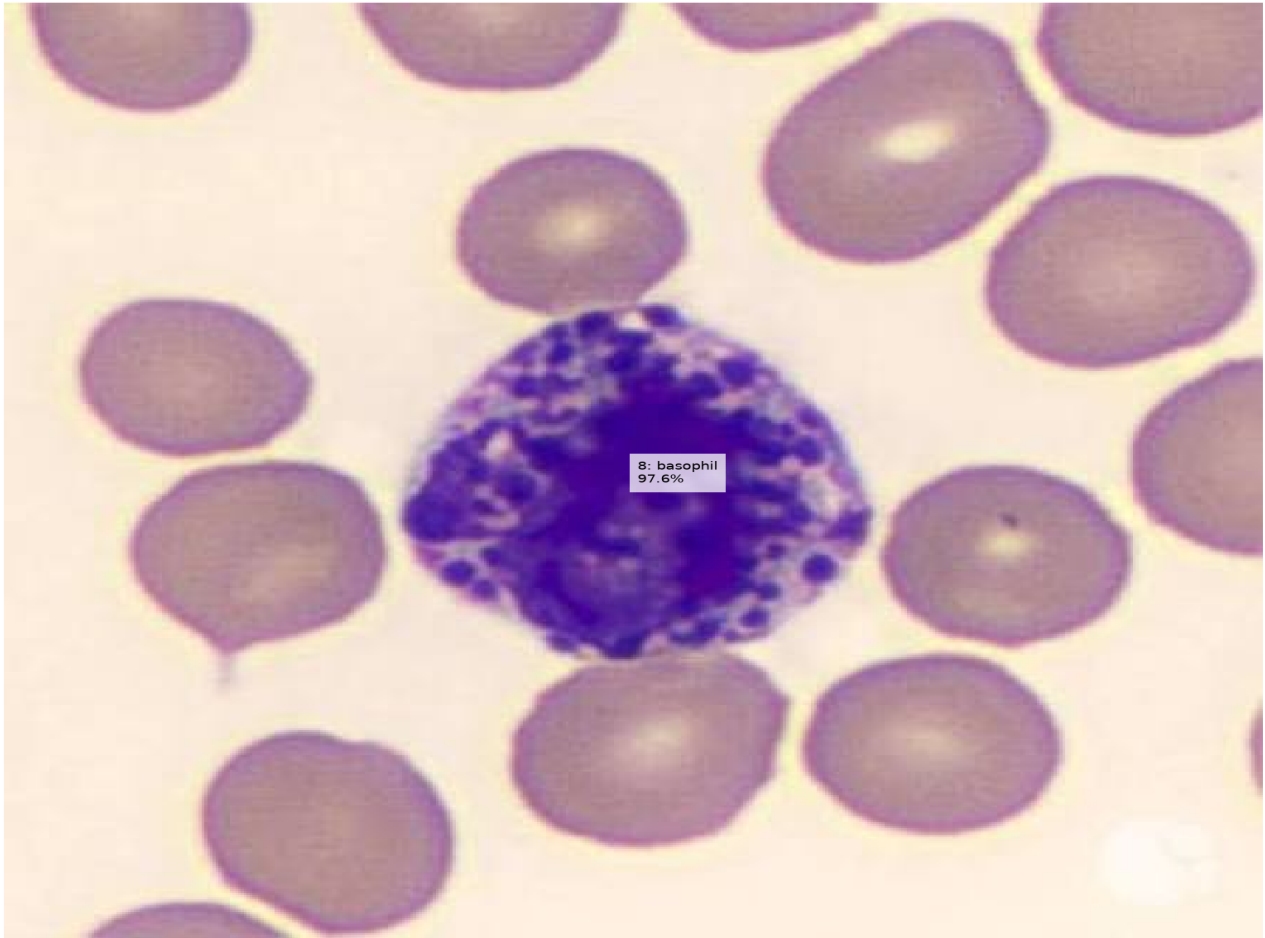

Cell-Level Results

ID	Area	Purple	Class	Confidence
8	20710	71.2	basophil	97.6%

Grad-CAM Explanations

Grad-CAM shows regions contributing to the selected model class. It is an explanation of model behavior, not a direct measurement of a biological structure.

SHAP Explanations

SHAP attribution magnitude summarizes the magnitude of pixel-level contribution toward the selected class. It is a model explanation, not a clinical measurement.

SHAP output was unavailable for this analysis.

Interpretation and Limitations

- The classifier has five closed-set classes: basophil, erythroblast, monocyte, myeloblast and segmented neutrophil.
- RBCs and unknown cell types are not explicit classes and can therefore be forced into one of the five outputs.
- Cellpose segmentation quality affects every downstream result.
- Stain, microscope, illumination and acquisition differences can cause domain shift.
- Confidence is not a calibrated probability of disease.
- A myeloblast-like prediction is not equivalent to an AML diagnosis.
- Clinical assessment requires appropriate hematopathology and laboratory testing.